

---

## Eye Disease Detection Using Deep Learning

### System Required:

---


Windows 8 machine Install with two web browser Bandwidth of 30mbps

### Project Description:

In this project we are classifying various types of Eye Diseases that people get due to various reasons like age, diabetes, etc. These diseases are majorly classified into 4 categories namely Normal, cataract, Diabetic Retinopathy & Glaucoma. Deep-learning (DL) methods in artificial intelligence (AI) play a dominant role as high-performance classifiers in the detection of the Eye Diseases using images.

Transfer learning has become one of the most common techniques that has achieved better performance in many areas, especially in image analysis and classification. We used Transfer Learning techniques like Inception V3, VGG19, exception V3 that are more widely used as a transfer learning method in image analysis and they are highly effective.

### Technical Architecture:

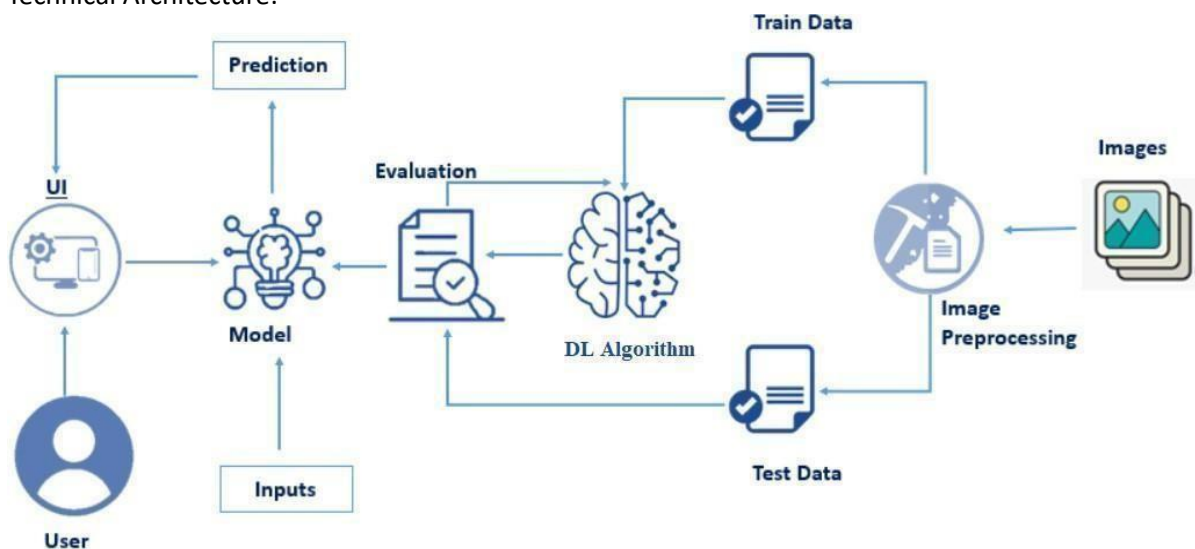

## Project Flow


- The user interacts with the UI (User Interface) to choose the image.
- The chosen image is analyzed by the model integrated with flask application.
- The Xception Model analyzes the image, then showcases the prediction on the Flask UI.

To accomplish this, we have to complete all the activities and tasks listed below

- Data Collection.
  - Create a Train and Test path.
- Data Pre-processing.
- Import the required library
- Configure ImageDataGenerator class
- Apply ImageDataGenerator functionality to Trainset and Testset
- Model Building
  - Pre-trained CNN model as a Feature Extractor
  - Adding Dense Layer
  - Configure the Learning Process
  - Train the model
  - Save the Model
  - Test the model
- Application Building
  - Create an HTML file

## Project Structure

Create a Project folder which contains files as shown below



Flask	19-01-2023 14:40
static	19-01-2023 14:17
templates	19-01-2023 15:06
uploads	19-01-2023 14:43
app.py	19-01-2023 14:40
evgg.h5	18-01-2023 10:55
IBM	18-01-2023 12:13
cnn_evgg.tgz	18-01-2023 12:11
eye_disease_detection_IBM.ipynb	18-01-2023 12:13
training	18-01-2023 12:01
evgg.h5	18-01-2023 10:55
eye_disease_detection.ipynb	18-01-2023 12:01
kaggle.json	17-01-2023 21:39
apikey.json	18-01-2023 12:02
Eye Disease Detection Using Deep Learning.docx	19-01-2023 16:12

- The Dataset folder contains the training and testing images for training our model.
- For building a Flask Application we need HTML pages stored in the **templates** folder, CSS for styling the pages stored in the static folder and a python script **app.py** for server side scripting
- The IBM folder consists of a trained model notebook on IBM Cloud.
- Training folder consists of eye\_disease\_detection.ipynb model training file & evgg.h5 is saved model

## Data Collection

There are many popular open sources for collecting the data. Eg: kaggle.com, UCI repository, etc.

### Download The Dataset

Collect images of Eye Diseases then organize into subdirectories based on their respective names as shown in the project structure. Create folders of types of Eye Diseases that need to be recognized.

In this project, we have collected images of 4 types of Eye Diseases images like Normal, cataract, Diabetic Retinopathy & Glaucoma and they are saved in the respective sub directories with their respective names.

You can download the dataset used in this project using the below link

Dataset:- <https://www.kaggle.com/datasets/gunavenkatdoddi/eye-diseases-classification>

**Note: For better accuracy train on more images**

We are going to build our training model on Google colab.

We will be connecting Kaggle with Google Colab because the dataset is too big to import using the following code:

```
!pip install -q kaggle

!mkdir ~/.kaggle

!cp kaggle.json ~/.kaggle/

!kaggle datasets download -d gunavenkatdoddi/eye-diseases-classification

Warning: Your Kaggle API key is readable by other users on this system! To fix this, you can run 'chmod 600 /root/.kaggle/kaggle.json'
Downloading eye-diseases-classification.zip to /content
 99% 736M/736M [00:07<00:00, 108MB/s]
100% 736M/736M [00:07<00:00, 101MB/s]
```

### Create Training And Testing Dataset

To build a DL model we have to split training and testing data into two separate folders. But in this project dataset folder training and testing folders are not present. So, in this case we have to separate the data into train & test folders.

```
import splitfolders

splitfolders.ratio('/content/dataset', output="output", seed=1337, ratio=(.8, 0.2))

Copying files: 4217 files [00:02, 1558.61 files/s]
```

Four different transfer learning models are used in our project and the best model (VGG19) is selected. The image input size of VGG19 model is 224, 224.

```
IMAGE_SIZE = [224, 224]
```

### Image Preprocessing

In this milestone we will be improving the image data that suppresses unwilling distortions or enhances some image features important for further processing, although perform some geometric transformations of images like rotation, scaling, translation, etc.

Link : [link](#)

### Importing The Libraries

Import the necessary libraries as shown in the image

```
import splitfolders
from tensorflow.keras.preprocessing.image import ImageDataGenerator
from PIL import ImageFile
ImageFile.LOAD_TRUNCATED_IMAGES = True
from tensorflow.keras.applications.vgg19 import VGG19, preprocess_input
from tensorflow.keras.layers import Flatten, Dense
from tensorflow.keras.models import Model
from tensorflow.keras.models import load_model
from tensorflow.keras.preprocessing import image
import numpy as np
import matplotlib.pyplot as plt
import matplotlib.image as mpimg
%matplotlib inline
```

## Configure ImageDataGenerator Class

ImageDataGenerator class is instantiated and the configuration for the types of data augmentation. There are five main types of data augmentation techniques for image data; specifically:

- Image shifts via the `width_shift_range` and `height_shift_range` arguments.
- The image flips via the `horizontal_flip` and `vertical_flip` arguments.
- Image rotations via the `rotation_range` argument.
- Image brightness via the `brightness_range` argument.
- Image zoom via the `zoom_range` argument.

An instance of the ImageDataGenerator class can be constructed for train and test.

```
from tensorflow.keras.preprocessing.image import ImageDataGenerator
train_datagen = ImageDataGenerator(rescale = 1./255,
                                   shear_range = 0.2, zoom_range = 0.2,
                                   horizontal_flip = True)

test_datagen = ImageDataGenerator(rescale = 1./255)
```

## Apply ImageDataGenerator Functionality To Train Set And Test Set

Let us apply ImageDataGenerator functionality to the Train set and Test set by using the following code. For Training set using `flow_from_directory` function.

This function will return batches of images from the subdirectories. Arguments:

- `directory`: Directory where the data is located. If labels are "inferred", it should contain subdirectories, each containing images for a class. Otherwise, the directory structure is ignored.
- `batch_size`: Size of the batches of data which is 64.
- `target_size`: Size to resize images after they are read from disk.
- `class_mode`:
  - 'int': means that the labels are encoded as integers (e.g. for `sparse_categorical_crossentropy` loss).
  - 'categorical' means that the labels are encoded as a categorical vector (e.g. for `categorical_crossentropy` loss).
  - 'binary' means that the labels (there can be only 2) are encoded as float32 scalars with values 0 or 1 (e.g. for `binary_crossentropy`).
  - None (no labels).

```
training_set = train_datagen.flow_from_directory(
    '/content/output/train',
                                target_size = (224, 224),
                                batch_size = 64,
                                class_mode = 'categorical')

test_set = test_datagen.flow_from_directory('/content/output/val',
                                target_size = (224, 224),
                                batch_size = 64,
                                class_mode = 'categorical')

Found 3372 images belonging to 4 classes.
Found 845 images belonging to 4 classes.
```

Total the dataset is having 3372 train images, 845 test images divided under 4 classes.

## Model Building

Now it's time to build our model. Let's use the pre-trained model which is VGG19, one of the convolution

neural net (CNN) architecture which is considered as a very good model for Image classification. Deep understanding on the VGG19 model – [Link](#) is referred to in the prior knowledge section. Kindly refer to it before starting the model building part.

### Pre-Trained CNN Model As A Feature Extractor

For one of the models, we will use it as a simple feature extractor by freezing all the five convolution blocks to make sure their weights don't get updated after each epoch as we train our own model.

Here, we have considered images of dimension (224,224,3).

Also, we have assigned `include_top = False` because we are using convolution layer for features extraction and wants to train fully connected layer for our image classification (since it is not the part of Imagenet dataset)

Flatten layer flattens the input. Does not affect the batch size.

```
VGG19 = VGG19(input_shape=IMAGE_SIZE + [3], weights='imagenet',include_top=False)

Downloading data from https://storage.googleapis.com/tensorflow/keras-applications/vgg19/vg
80134624/80134624 [=====] - 1s 0us/step

for layer in VGG19.layers:
    layer.trainable = False

x = Flatten()(VGG19.output)
```

## Adding Dense Layers

```
x = Flatten()(VGG19.output)

prediction = Dense(4, activation='softmax')(x)

model = Model(inputs=VGG19.input, outputs=prediction)
```

A dense layer is a deeply connected neural network layer. It is the most common and frequently used layer. Let us create a model object named model with inputs as VGG19.input and output as dense layer.

The number of neurons in the Dense layer is the same as the number of classes in the training set. The neurons in the last Dense layer, use softmax activation to convert their outputs into respective probabilities. Understanding the model is a very important phase to properly use it for training and prediction purposes.

Keras provides a simple method, summary to get the full information about the model and its layers.

## Configure The Learning Process

The compilation is the final step in creating a model. Once the compilation is done, we can move on to the training phase. The loss function is used to find errors or deviations in the learning process. Keras requires a loss function during the model compilation process.

Optimization is an important process that optimizes the input weights by comparing the prediction and the loss function. Here we are using adam optimizer



Metrics are used to evaluate the performance of your model. It is similar to the loss function, but not used in the training process

```
model.compile(loss='categorical_crossentropy', optimizer='adam', metrics=['accuracy'])
```

## Train The Model

Now, let us train our model with our image dataset. The model is trained for 50 epochs and after every epoch, the current model state is saved if the model has the least loss encountered till that time. We can see that the training loss decreases in almost every epoch and probably there is further scope to improve the model.

**.fit** functions used to train a deep learning neural network

### Arguments:

- **steps\_per\_epoch**: it specifies the total number of steps taken from the generator as soon as one epoch is finished and the next epoch has started. We can calculate the value of **steps\_per\_epoch** as the total number of samples in your dataset divided by the batch size.
- **Epochs**: an integer and number of epochs we want to train our model for.
- **validation\_data** can be either:
  - an inputs and targets list
  - a generator
  - an inputs, targets, and sample\_weights list which can be used to evaluate the loss and metrics for any model after any epoch has ended.
- **validation\_steps**: only if the **validation\_data** is a generator then only this argument can be used. It specifies the total number of steps taken from the generator before it is stopped at every epoch and its value is calculated as the total number of validation data points in your dataset divided by the validation batch size.

```
r = model.fit(
    training_set,
    validation_data=test_set,
    epochs=50,
    steps_per_epoch=len(training_set),
    validation_steps=len(test_set)
)
```

```
Epoch 33/50
53/53 [=====] - 64s 1s/step - loss: 0.2517 - accuracy: 0.9015 - val_loss: 0.3413 - val_accuracy: 0.8817
Epoch 34/50
53/53 [=====] - 64s 1s/step - loss: 0.2747 - accuracy: 0.8986 - val_loss: 0.6687 - val_accuracy: 0.7787
Epoch 35/50
53/53 [=====] - 65s 1s/step - loss: 0.2640 - accuracy: 0.8923 - val_loss: 0.5145 - val_accuracy: 0.8154
Epoch 36/50
53/53 [=====] - 65s 1s/step - loss: 0.2475 - accuracy: 0.9036 - val_loss: 0.6416 - val_accuracy: 0.7893
Epoch 37/50
53/53 [=====] - 64s 1s/step - loss: 0.2497 - accuracy: 0.9021 - val_loss: 0.3830 - val_accuracy: 0.8781
Epoch 38/50
53/53 [=====] - 64s 1s/step - loss: 0.2462 - accuracy: 0.9054 - val_loss: 0.5982 - val_accuracy: 0.8118
Epoch 39/50
53/53 [=====] - 65s 1s/step - loss: 0.2618 - accuracy: 0.8980 - val_loss: 0.4216 - val_accuracy: 0.8710
Epoch 40/50
53/53 [=====] - 64s 1s/step - loss: 0.2343 - accuracy: 0.9012 - val_loss: 0.5007 - val_accuracy: 0.8213
Epoch 41/50
53/53 [=====] - 64s 1s/step - loss: 0.2198 - accuracy: 0.9075 - val_loss: 0.4164 - val_accuracy: 0.8580
Epoch 42/50
53/53 [=====] - 64s 1s/step - loss: 0.2489 - accuracy: 0.9024 - val_loss: 0.6112 - val_accuracy: 0.8107
Epoch 43/50
53/53 [=====] - 64s 1s/step - loss: 0.2450 - accuracy: 0.9021 - val_loss: 0.3405 - val_accuracy: 0.8888
Epoch 44/50
53/53 [=====] - 65s 1s/step - loss: 0.2369 - accuracy: 0.9101 - val_loss: 0.4713 - val_accuracy: 0.8343
Epoch 45/50
53/53 [=====] - 64s 1s/step - loss: 0.2593 - accuracy: 0.8983 - val_loss: 0.3722 - val_accuracy: 0.8722
Epoch 46/50
53/53 [=====] - 64s 1s/step - loss: 0.2262 - accuracy: 0.9078 - val_loss: 0.4198 - val_accuracy: 0.8615
Epoch 47/50
53/53 [=====] - 64s 1s/step - loss: 0.2300 - accuracy: 0.9090 - val_loss: 0.4275 - val_accuracy: 0.8639
Epoch 48/50
53/53 [=====] - 65s 1s/step - loss: 0.2230 - accuracy: 0.9137 - val_loss: 0.3193 - val_accuracy: 0.8959
Epoch 49/50
53/53 [=====] - 64s 1s/step - loss: 0.2204 - accuracy: 0.9081 - val_loss: 0.3112 - val_accuracy: 0.8970
Epoch 50/50
53/53 [=====] - 64s 1s/step - loss: 0.2597 - accuracy: 0.8956 - val_loss: 0.4516 - val_accuracy: 0.8355
```

From the above run time, we can observe that at 50<sup>th</sup> epoch the model is giving the best accuracy.

### Testing the model:

Evaluation is a process during the development of the model to check whether the model is the best fit for the given problem and corresponding data. Load the saved model using `load_model`.

```
from tensorflow.keras.models import load_model
from tensorflow.keras.preprocessing import image
import numpy as np
model = load_model("/content/evgg.h5")
```

Taking an image as input and checking the results

```

img = image.load_img(r"/content/output/val/normal/2365_left.jpg",target_size= (224,224))#loading of the image
x = image.img_to_array(img)#image to array
x = np.expand_dims(x,axis = 0)#changing the shape
preds=model.predict(x)
pred=np.argmax(preds,axis=1)
index=['cataract','diabetic_retinopathy','glaucoma','normal']
result=str(index[pred[0]])
result

1/1 [=====] - 1s 806ms/step
'normal'

```

So our model has predicted the label correctly as Nomal.

## Application Building

In this section, we will be building a web application that is integrated to the model we built. A UI is provided for the uses where he has to enter the values for predictions. The enter values are given to the saved model and prediction is showcased on the UI.

This section has the following tasks

- Building HTML Pages
- Building server side script

## Building Html Pages

For this project create one HTML file namely

- index.html


Let's see how our index.html page looks like:

predict section

We're **determined** for  
your **better life.**

You can get the care you need 24/7 – be it online or in person. You will be treated by caring specialist doctors.

## Make an Appointment

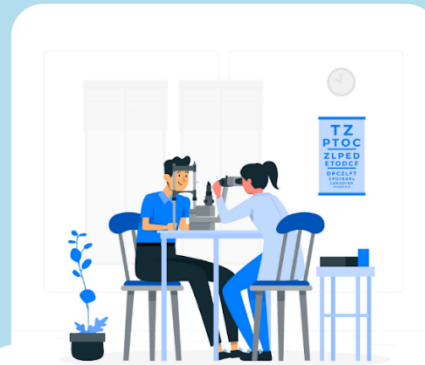

**Eye Care with Top Professionals  
and In Budget.**

We've built a healthcare system that puts your needs first. For us, there is nothing more important than the health of you and your loved ones.

[Learn more](#)


## ABOUT US

**We are developing a healthcare system around you**

We think that everyone should have easy access to excellent healthcare. Our aim is to make the procedure as simple as possible for our patients and to offer treatment no matter where they are — in person or at their convenience.

[Learn more](#)



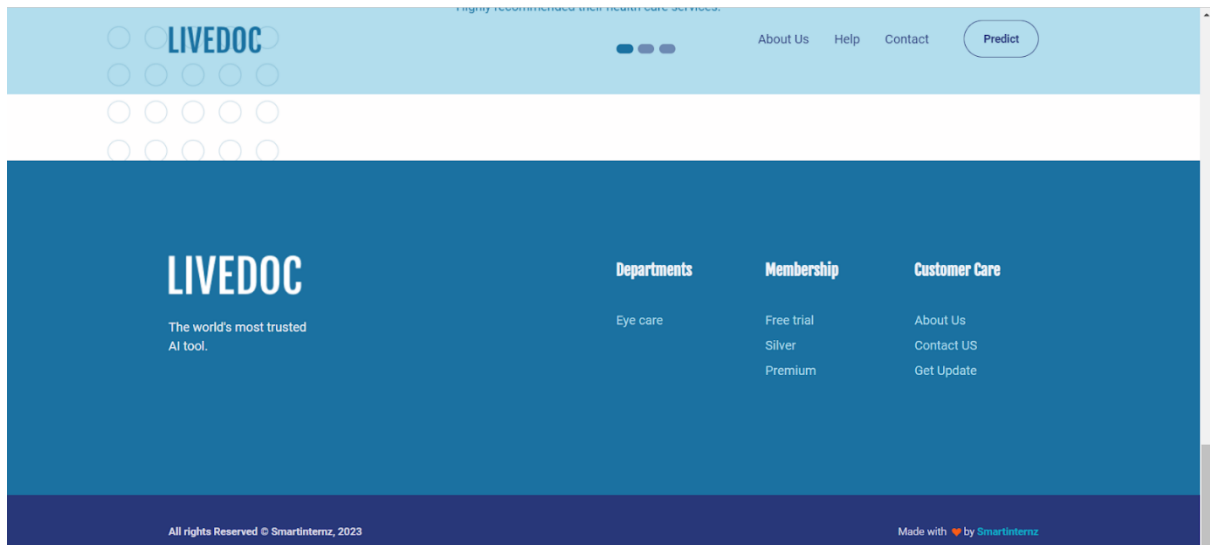

- Eye Disease Detection Using Deep Learning
  - Prerequisites
  - Prior Knowledge
  - Project Objectives
  - Project Flow
  - Project Structure
  - Data Collection
    - Download The Dataset
    - Create Training And Testing Dataset
    - Image Preprocessing
  - Image Preprocessing
    - Importing The Libraries
    - Configure ImageDataGenerator Class
    - Apply ImageDataGenerator Functionality To Train Set And Test Set
  - Model Building
    - Pre-Trained CNN Model As A Feature Extractor
    - Adding Dense Layers
    - Configure The Learning Process
    - Train The Model
  - Save The Model
  - Application Building
    - Building Html Pages
    - Build Python Code
    - Run The Application
  - Output

### **Build Python Code**

Import the libraries

```
import numpy as np
import os
from flask import Flask, request, render_template
from tensorflow.keras.models import load_model
from tensorflow.keras.preprocessing import image
from tensorflow.keras.applications.resnet_v2 import preprocess_input
```

Loading the saved model and initializing the flask app

```
model=load_model(r"evgg.h5")

app=Flask(__name__)
```

Render HTML pages:

```
@app.route('/')
def index():
    return render_template('index.html')

@app.route('/home')
def home():
    return render_template("index.html")

@app.route('/inp')
def inp():
    return render_template("img_input.html")
```

Once we uploaded the file into the app, then verifying the file uploaded properly or not. Here we will be using declared constructor to route to the HTML page which we have created earlier.

In the above example, '/' URL is bound with index.html function. Hence, when the home page of the web server is opened in browser, the html page will be rendered. Whenever you enter the values from the html page the values can be retrieved using POST Method.

```

@app.route('/predict', methods=["GET", "POST"])
def res():
    if request.method=="POST":
        f=request.files['image']
        basepath=os.path.dirname(__file__)
        filepath=os.path.join(basepath, 'uploads', f.filename)
        f.save(filepath)

        img=image.load_img(filepath, target_size=(224, 224, 3))
        x=image.img_to_array(img)
        x=np.expand_dims(x, axis=0)

        img_data=preprocess_input(x)
        prediction=np.argmax(model.predict(img_data), axis=1)

        index=['cataract', 'diabetic_retinopathy', 'glaucoma', 'normal']

        result=str(index[prediction[0]])
        print(result)
        return render_template('output.html', prediction=result)

```

Here we are routing our app to predict function. This function retrieves all the values from the HTML page using Post request. That is stored in an array. This array is passed to the model.predict() function. This function returns the prediction. And this prediction value will be rendered to the text that we have mentioned in the index.html page earlier.

#### Main Function:

```

if __name__ == "__main__":
    app.run(debug=False)

```

#### Run The Application

- Open Spyder
- Navigate to the folder where your Python script is.
- Now click on the green play button above.
- Click on the predict button from the top right corner, enter the inputs, click on the Classify button, and see the result/prediction on the web.




```
* Serving Flask app "app" (lazy loading)
* Environment: production
  WARNING: This is a development server. Do not use it in a
production deployment.
  Use a production WSGI server instead.
* Debug mode: off
* Running on http://127.0.0.1:5000/ (Press CTRL+C to quit)
```

## Output

The home page looks like this. When you click on the Predict button, you'll be redirected to the predict section

---


### Prediction Result

 Uploaded Image

**Predicted Label: normal**

Upload Another Image

Input :

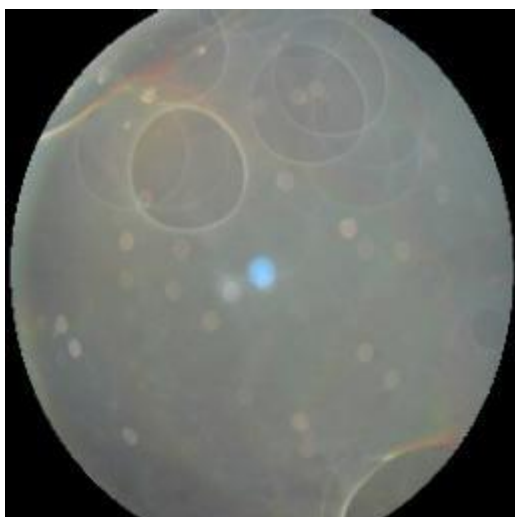

Once you upload the image and click on Predict button, the output will be displayed in the below page

Output :

### Prediction Result

 Uploaded Image

**Predicted Label: glaucoma**

Upload Another Image